

## THE 500KLHB

# PRIMARY CARE DIGITAL SCALE WITH LIVE HANDLEBARS

The number-one-selling physician scale in North America, the all-digital Health o meter® 500KLHB primary care scale provides the optimum in quality, reliability, and ease of use. The unique live handlebars provide stability that many patients need, without affecting the accuracy of the weight reading. Designed to meet the needs of doctor's offices, hospitals, clinics, and acute and long-term care facilities, the scale features a built-in height rod in the scale column, a large non-skid platform, and wheels that enable portability. The scale gives providers and patients an eye-level weight indicator for patients up to 550 pounds, with a user-friendly digital interface that includes BMI, hold/release, zero, tare, reweigh, recall, KG/LB, clear, and enter functions. The 500KL offers the best in its class extra-wide 8" wheelbase making transport easy. The exclusive Everlock® function allows the measuring unit (LB or KG) to be permanently locked on the scale.

The first healthcare scale to provide easy integration and seamless communication with today's EMR applications, the 500KLHB includes the option to add wireless connectivity to interface with a computer, EMR software, or other electronic device.

With its low platform height, the 500KLHB is easy to set up and offers patients a more comfortable and stable weighing experience.

Product Code:

HOM500KLHB- Digital Scale with Handlebars

## THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers on our scale, it's not surprising that the 500KLHB is the most widely used scale in physician offices nationwide.


**150,000+**  
INSTALLED  
NATIONWIDE

**120+** LOCAL  
FIELD SPECIALISTS

**8"** WIDE WHEELBASE  
BEST IN ITS CLASS

**#1** TRUSTED  
PHYSICIAN SCALE

**100** YEARS  
OF SERVICE

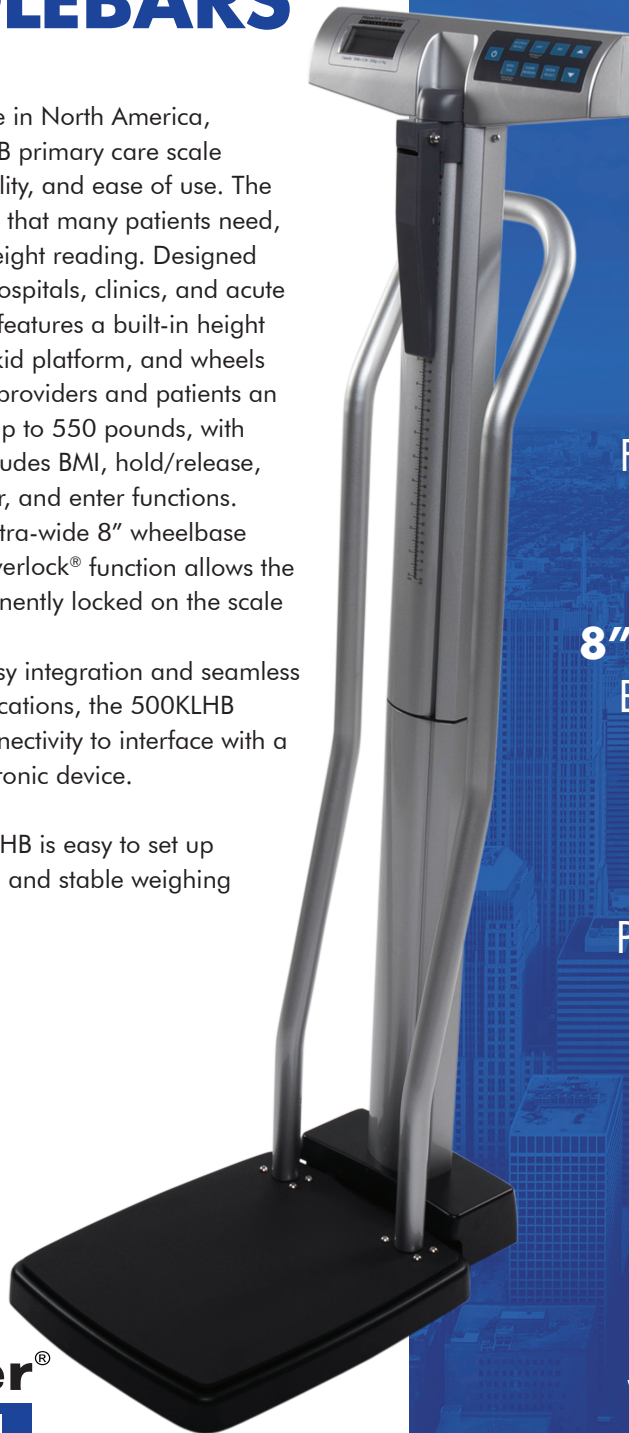

**Health o meter®**  
Professional

Health o meter® Professional 500KLHB

weigh**easier**®

## TECHNICAL SPECIFICATIONS

### Capacity

550 lb / 250 kg

### Resolution

0.2 lb / 0.1 kg

### EMR connectivity

Via optional Pelstar®  
wireless technology

### Display

1" / 25 mm, LCD display

### Batteries

6 AA batteries included,  
optional adapter  
(order# ADPT31)

### Functions

LB/KG conversion  
LB/KG lock out  
BMI  
Zero  
Tare, Pre-tare  
Hold/Release  
Auto Zero  
Auto Off  
Reweigh  
Recall  
Parent/Child weighing function  
Everlock® - permanent unit lock

### Platform Size (w x d x h)

13 ¾" x 16 ½" x 2 ⅜"  
(349 mm x 419 mm x 60 mm)  
Dent and rust proof platform

### Product Footprint

14 ⅜" x 21 ⅛" x 53 ⅝"  
(359 mm x 537 mm x 1362 mm)

### Product Weight

31 lb / 15 kg

### Height Rod

30"-84 ¼" / 76 cm-214 cm  
Graduation: ⅛" / 1 mm

### Portability

2 wheels, wide 8" wheelbase for ease  
of movement

### Warranty

2-year limited manufacturer warranty

## ORDERING INFORMATION

MODEL	DESCRIPTION
500KLHB	500KL with handlebar accessory
SS-500KL	ScaleSurance® 2 year extended warranty for 500 series*
ADPT31	Optional adapter for 500KLHB
PLUGSET30	Plug adapters for use outside North America

The handlebar is easy to attach to any new scale using the pre-drilled mounting holes and can also be attached to any existing scale using the easy-to-use drill guides.



**Health o meter®**  
Professional



\*Not available in all countries.

500KLHB\_F&B\_20220921