

Diagnostic Instruments



ATLANTIC
MEDICAL
SOLUTIONS



Our line of diagnostic EENT instruments totals nearly 900 model options, the largest selection by unique part number in the United States. We offer two brands: economical 2.5V portable Proscope™ models, and feature-rich 2.5V and 3.5V Diagnostix™ instruments in pocket, portable, desk, and wall systems.

Choose from our line of full-size 3.5V and 2.5V diagnostic instruments or compact pocket sets, depending on clinical needs or budget. We offer four different instrument heads, for virtually every type of physical exam:

- **Otoscopes for examination of the inner ear**
- **Ophthalmoscopes for examination of the eye**
- **Throat illuminators for examination of the mouth and throat**
- **Dermoscopes for examination of the skin.**
- **Our 2.5V otoscope is available in a veterinary model, as well.**



PMV
Otoscope
Page 5



Standard
Otoscope
Page 5



Dermoscope
Page 7



Ophthalmoscope
Page 6



Throat
Illuminator
Page 7

Purchase instruments individually or choose one of our popular combination sets, with hard and soft cases for storage.



ADC is one of a few firms to hold our own 510(k) and be classed as a manufacturer of EENT instruments

All ADC otoscopes accommodate Welch Allyn® KleenSpec® disposable ear specula for the ultimate in convenience and affordability.



SOURCED GLOBALLY. INSPECTED, ASSEMBLED and PACKAGED LOCALLY

We contract with the world's leading manufacturers in Pakistan, Germany, and the U.S.A. to produce components to our exacting specifications. Then, using custom-built inspection equipment and factory-trained personnel, we perform rigorous testing in our New York headquarters before the instruments are assembled and packaged.



COMPREHENSIVE WARRANTIES

As with all ADC warranties, coverage for our Diagnostix and Proscope brand EENT instruments is among the industry's most comprehensive.

- **Proscope 2.5V Portable instruments are warranted for one year**
- **Diagnostix 2.5V Pocket instruments are warranted for two years**
- **Diagnostix 3.5V instrument heads, handles, and desk chargers are warranted for two years**
- **Diagnostix Adstation wall transformers and extension modules are warranted for five years**
- **Optical components are warranted for life**
- **Diagnostix AdLED lamps are warranted for life**

Special Medical Student Warranty: Diagnostix 3.5V Portable sets are covered through the second year of residency.

Power Options for Any Clinical Need

Our 2.5V portable and pocket instruments are available with alkaline battery handles in pocket (AA) and portable (C) sizes.

Our 3.5V rechargeable battery handles (portable and desktop platforms) utilize Li-ion battery technology – the same found in smart phones, tablets, and electric cars. Combined with super-efficient LED lamps, our Li-ion battery provides up to 50 hours continuous operation between charges (compared to just two to three for NiMH or NiCad) and a whopping 1,200 charging cycles – 2.5X that of NiMH or NiCad. That means nearly four years of daily recharging, even longer if charged less often.

ADC's 3.5V Portable Diagnostix line has a **new USB rechargeable handle option!**
See page 8



Build Your Own

Our online Diagnostix Builder lets you create your perfect portable, desk, or wall configuration. Compare options and mix and match the features you're looking for.

www.adctoday.com/diagnostix/build



ADC Adstations include a universal power supply with US plug. For international use, our Global Adapter Kit (item 5660TN-GAK) can be ordered separately. The kit includes three replacement plugs for the EU, UK, and Australia that connect to our universal power supply.



Diagnostic Adstation
See pages 14-15

A Range of Illumination Options

Our **portable 2.5V instruments** are available with low-cost vacuum lamps.

Pocket instruments offer xenon illumination for a brighter, whiter light (3,200K), or AdLED lighting for true tissue color (4,200K).

Our **3.5V instruments** are available with either xenon lamps or our third generation LED technology — known as AdLED — for the ultimate in performance.



Adtronic™ Light Control


Our 2.5V portable instruments offer mechanical rheostatic light control. Our 3.5V instruments feature patented **Adtronic Electronic Rheostatic** light output control that provides infinite, continuously variable control. And because it's electronic, we built in an auto-off feature that engages after two minutes of inactivity to further conserve battery strength and lamp life.



Pocket Sets: Available in six colors!

ADC's Pocket Series Diagnostic Instruments include both an Otoscope and Ophthalmoscope. Choice of LED or xenon illumination and a compact size make these indispensable instruments for today's busy practitioner.

Comparing 3.5V Lighting Options

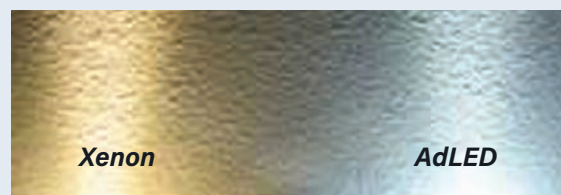

Xenon provides a bright white light (3,200K) with cooler operating temperatures.

AdLED provides near-perfect performance with a color temperature of 4,200K, 17,000 lux, 90+ color rendering index (CRI), and a life expectancy of 100,000 hours (compared with 20-30 hours for xenon). When combined with our Li-ion battery technology AdLED bulbs extend performance to levels never before achieved in an EENT instrument.

Diagnostix AdLED lights are warranted for life.



Diagnostix 5110E:
Our compact, single-handle
pocket set.
See page 12

Instrument Heads

OTOSCOPES

Diagnostix™ 3.5V PMV (Panoramic Magnified View) Otoscope Head - Item 5422



- Twice the field of view and 66% greater magnification (11 mm field of view, 5X magnification) compared to standard otoscopes
- Significantly improved view of the auditory canal and tympanic membrane
- Fiber-optic light transmission for an unobstructed view
- Choose xenon or high-performance AdLED illumination
- Reflection-optimized, high-performance glass optics
- Adjustable focal depth
- Suitable for throat illumination and pneumatic otoscopy (insufflator sold separately)
- Compatible with Welch Allyn specula
- Anti-theft locking feature is standard
- Works with ADC's 5460, 5463, 5560, and 5660 handles (universal bayonet mount also compatible with other brands of 3.5V handles)
- Two-year warranty on instrument head
- Lifetime warranty on LED lamp
- Manufactured in Germany, inspected and packaged in the U.S.A.
- Specify "X" for xenon or "L" for LED when ordering

Diagnostix™ 3.5V Otoscope Instrument Head - Item 5420



- Fiberoptic light transmission for an unobstructed view
- Optical glass reduces reflection and provides 3x magnification
- Bidirectional swivel lens allows instrumentation
- Insufflator fitting for pneumatic otoscopy
- Compatible with Welch Allyn specula
- Instrument head uses universal bayonet mount for compatibility with other brands of 3.5V handles
- Anti-theft locking feature is standard
- Works with 5460, 5463, 5560, and 5660 handles
- Choose xenon or high-performance AdLED illumination
- Manufactured in Germany, inspected and packaged in the U.S.A.
- Two-year warranty on instrument
- Lifetime warranty on LED lamp
- Specify "X" for xenon or "L" for LED when ordering

Diagnostix™ 2.5V Pocket Otoscope Instrument Head - Item 5120N



- Choose AdLED or 2.5V xenon lamp for true tissue color
- Fiber optic light transmission for an unobstructed view
- Wide-angle swivel 2.5X viewing lens permits instrumentation
- Insufflator port for pneumatic otoscopy (fitting and insufflator sold separately)
- Works with 5160N handle
- Two-year warranty on instrument
- Manufactured in Pakistan, inspected and packaged in the U.S.A.
- Available in black, green (GR), royal blue (RB), purple (V), burgundy (BD), or gray (G)
- Specify "L" for LED when ordering

Proscope™ 2.5V Otoscope Instrument Head - Item 5220



- 2.5V vacuum lamp
- Direct illumination
- Bayonet locking mount
- Removable 3X magnifying lens
- Includes autoclavable, polypropylene ear specula in sizes 2.5mm, 3.5mm, and 4.5mm
- Welch Allyn disposable specula adapter included
- Works with 5211-5 battery handle
- One-year warranty on instrument
- Manufactured in Pakistan, inspected and packaged in the U.S.A.

OPHTHALMOSCOPES




Diagnostic™ Coax and Coax Plus 3.5V Ophthalmoscope Instrument Heads

- High-performance aspherical condenser lens reduces corneal reflection
- Ergonomic design and integral thumb rest improves stability and allows one-finger control of all functions
- Dust-proof housing and soft synthetic rubber brow rest for use with eyeglasses
- Magnified/illuminated aperture viewing window
- Works with ADC's 5460, 5463, 5560, and 5660 handles; universal bayonet mount also compatible with other brands of 3.5V handles
- Anti-theft locking feature is standard
- Choose xenon or high-performance AdLED
- Two-year warranty on instrument
- Lifetime warranty on LED lamp
- Specify "X" for xenon or "L" for LED when ordering

- Twenty-four aperture/filter combinations for improved versatility

Coax Plus Ophthalmoscope (Item 5442) also includes:

- Eighty-nine viewing lenses
Plus: single step from 1-45; Minus: single step from 1-44
- Secondary focusing wheel with zero retention for quick correction of ametropias, provides virtually unlimited fine focus or fast focus in lieu of the 89 viewing lenses
- Twenty-eight aperture/filter combinations for improved versatility

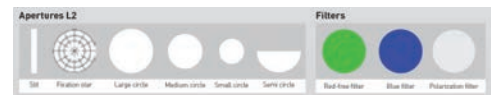

5440



Coax Ophthalmoscope (Item 5440) also includes:

- Twenty-nine viewing lenses from -35 to +40
Plus: 1-10, 12, 15, 20, 40; Minus: 1-10, 15, 20, 25, 30, 35



5442



Diagnostic™ 2.5V Pocket Ophthalmoscope Instrument Head - Item 5140N

- Choose AdLED or 2.5V xenon lamp for true tissue color
- Five aperture selections: small spot, large spot, semicircle, fixed, and red-free filter
- Nineteen viewing lenses from -20 to +20 diopters for better resolution
- Illuminated magnified lens dial
- Soft synthetic rubber brow rest is helpful when using with eyeglasses

- Meets or exceeds ISO 10942, ISO 15004-1, and ISO 15004-2
- Works with 5160N handle
- Two-year warranty on instrument
- Manufactured in Pakistan, inspected and packaged in the U.S.A.
- Available in black, green (GR), royal blue (RB), purple (V), burgundy (BD), or gray (G)
- Specify "L" for LED when ordering



Proscope™ 2.5V Ophthalmoscope Instrument Head - Item 5240

- 2.5V vacuum lamp
- 24 corrective lenses from -25 to +40 diopters
- Medium spot aperture
- Bayonet locking mount

- Works with 5211-5 handle
- One-year warranty
- Manufactured in Pakistan, inspected and packaged in the U.S.A.

THROAT ILLUMINATORS



Diagnostix™ 3.5V Throat Illuminator Instrument Head - Item 5426

- PMMA light transmission
- Tongue blade holder facilitates examination of the mouth and throat
- Manufactured in Germany, inspected and packaged in the U.S.A.
- Eject button discards used tongue blade
- Works with 5460, 5560, and 5660 handles
- Two-year warranty on instrument
- Lifetime warranty on LED lamp

Specify "X" for xenon or "L" for LED when ordering

DERMASCOPES



Diagnostix™ 3.5V Dermoscope Instrument Head - Item 5450

- Direct Illumination
- Optical grade 10X focusable magnifying lens
- Large knurled focusing ring
- Removable contact plate with mm scale
- Soft synthetic rubber eyepiece is helpful when using with eyeglasses
- Anti-theft locking feature is standard
- LED illumination only
- Two-year warranty on instrument
- Lifetime warranty on LED lamp



Proscope™ 2.5V Dermoscope Instrument Head - Item 5312-3

- Bright 2.5V halogen illumination ensures true tissue color for accurate identification
- 10X optic system
- Removable contact plate is autoclavable
- Large knurled focusing ring
- Reticle contact plate etched onto the lens in millimeters permits easy measurement of suspicious lesions
- Bayonet locking mount
- Soft synthetic rubber eyepiece is helpful when using with eyeglasses
- Works with 5260 handle
- One-year warranty on instrument
- Manufactured in Pakistan, inspected and packaged in the U.S.A.

Power Supply



Global Adapter Kit

Diagnostix™ 3.5V Global Adapter Kit

Our Global Adapter Kit is a set of three AC plugs that can convert our Diagnostix 3.5V universal power supply for use outside the U.S.

- Includes three replacement plugs for the EU, UK, and Australia (Power cord not included.)
- Plugs quickly connect to the universal power supply that ships with Diagnostix Adstation wall transformers and extension modules, as well as the Diagnostix desk set charging station
- Two-year warranty

Ordering Information

Item	Description
5660TN-GAK	Global Adapter Kit, EU, UK, AU



Diagnostix™ 3.5V Plug-In Rechargeable Handle

- C-sized handle with rechargeable 1400mAh Li-ion battery
- Top half of handle unscrews to access wall plug (US and Euro fittings available)
- Li-ion battery provide 2.5 hours continuous use with Xenon, 6.5 hours use with LED when charged for 24 hours
- Bayonet mount works with all ADC 3.5V heads (5420, 5442, 5440, 5442, 5450, and 5426) and is compatible with all other major brands
- Handle is CNC machined from solid brass then chrome plated in black with a ribbed surface for a sure grip
- Adtronic™ electronic rheostatic control provides infinite light control and auto-off
- Two-year warranty



Ordering Information

Item	Description
5460BAT	Replacement battery for 5460



Ordering Information

Item	Description
5463BAT	Replacement battery for 5463



Diagnostix™ 3.5V USB Rechargeable Handle

- C sized handle with USB rechargeable 18650 Li-ion battery, 3.5V, 2600 mAh, 9.62Wh
- USB-C charging cord included
- Li-ion battery provide 3.5 hours continuous use with Xenon, 10 hours use with LED when charged for 24 hours
- Bayonet mount works with all ADC 3.5V heads (5420, 5442, 5440, 5442, 5450, and 5426) and is compatible with all other major brands
- Handle is CNC machined from solid brass then chrome plated in black with a ribbed surface for a sure grip
- Adtronic™ electronic rheostatic control provides infinite light control and auto-off
- Two-year warranty

Diagnostix™ 3.5V Desk Rechargeable Handle

- C-sized handle with rechargeable 2600mAh Li-ion battery
- Requires Desk Charger Base, Item 5500, sold separately
- Li-ion battery provides 2-3 hours of continuous use with xenon lamps and up to 50 hours of continuous use with LED lamps
- Bayonet mount works with all ADC 3.5V heads (5420, 5442, 5440, 5442, 5450, and 5426) and is compatible with all other major brands
- Handle is CNC machined from solid brass then chrome plated in black with a ribbed surface for a sure grip
- Adtronic™ electronic rheostatic control provides infinite light control and auto-off
- Two-year warranty



Ordering Information

Item	Description
5560BAT	Replacement battery for 5560



Desk Recharger

Ordering Information

Item	Description
5500	3.5V Diagnostix™ Desk Recharger
5660TN-GAK	Global Adapter Kit, EU, UK, AU



Diagnostix™ 3.5V Desk Recharger

- Compact and portable desk charging station for two Diagnostix™ rechargeable 5560 handles
- Includes reducing sleeve for use with smaller AA sized handles
- Colored LED power indicator displays charging state
- Storage compartment accommodates a small selection of specula or spare lamps
- Includes wall-mount bracket for permanently mounting to wall or cabinet
- Includes universal power supply with U.S. plug; for international use, our Global Adapter Kit (5660TN-GAK) can be ordered separately



Diagnostix™ 2.5V Pocket Battery Handle

- AA 2.5V battery handle
- On/Off power switch integrated into pocket clip
- Metal-lined thermoplastic barrel
- Two-year warranty
- Manufactured in Pakistan, inspected and packaged in the USA
- Works with 5120N and 5140N instrument heads
- Available in black, green (GR), royal blue (RB), purple (V), burgundy (BD), or gray (G)

Ordering Information

Item	Description
5160N*	Diagnostix™ 2.5V Pocket Battery Handle



* Specify "BD" for Burgundy, "G" for Gray, "GR" for Green, "RB" for Royal Blue, or "V" for Purple when ordering. (Black is standard.)

Proscope™ 2.5V Portable Handles

- C-cell 2.5V power supply
- Rheostatic illumination control
- Knurled brass handle for positive grip
- Bayonet locking mount
- One-year warranty
- Manufactured in Pakistan, inspected and packaged in the USA
- 5211-5 works with 5220, 5240, and 5215-TH instrument heads
- 5260 works with 5312-3 instrument head

Ordering Information

Item	Description
5211-5	Proscope™ 2.5V Portable Handle / Oto, Oph, Throat Illum.
5260	Proscope™ 2.5V Portable Handle / Dermoscope



5660T



5660E

Diagnostix™ Adstation™ Wall Transformer and Wall Extension Module

- Wall-mount power supply helps improve workflow
- Start with a single power handle and add up to three extension modules
- Magnetic switch activates power when handle is removed from cradle, automatically powers off when handle is replaced
- Adtronic™ electronic rheostatic light-intensity control provides infinite continuous control
- Adtronic™ auto-off feature after two minutes provides additional safety if instrument is mis-cradled
- Easy to install with included bracket or may be mounted to our wallboard systems
- Bayonet instrument mount is compatible with all major brands of 3.5V instrument heads
- Includes universal power supply with U.S. plug. For international use, our Global Adapter Kit (5600TN-GAK) can be ordered separately
- Five-year warranty
- Works with 5420, 5422, 5440, 5442, 5450, and 5426 instrument heads

Ordering Information

Item	Description
5660T	Diagnostix™ Adstation™ Wall Transformer
5660E	Diagnostix™ Adstation™ Wall Extension Module
5660TN-GAK	Global Adapter Kit, EU, UK, AU





Diagnostix™ Desk Sets

The Diagnostix™ Desk Set features

- Compact and portable desk charging station for two Diagnostix™ rechargeable handles (5560)
- Reducing sleeve to allow use with smaller AA-sized handles
- Colored LEDs indicate power and charging state
- Handles charge fully overnight, providing three hours continuous use with xenon lamps or 50 hours continuous use with LED lamps
- Storage compartment accommodates a small selection of specula or spare lamps
- Includes wall-mount bracket for permanently mounting to wall or cabinet
- Includes universal power supply with U.S. plug. For international use, our Global Adapter Kit (5660TN-GAK) can be ordered separately

Diagnostix™ Desk Set Model Descriptions



Component Model	5422 PMV Otoscope	5420 Standard Otoscope	5440 Coax Ophthal.	5442 Coax Plus Ophthal.	5560 Handle (x2)	5500 Desk Charger
5580*	✓		✓		✓	✓
55802*	✓			✓	✓	✓
5510*		✓	✓		✓	✓
55102*		✓		✓	✓	✓

* Indicates lamps: X = Xenon, L = LED, A = LED Otoscope & Xenon Ophthalmoscope, B = Xenon Otoscope & LED Ophthalmoscope

Ordering Information

Item	Description
5510*	Diagnostix™ Desk Set, Otoscope, Coax Ophthal
55102*	Diagnostix™ Desk Set, Otoscope, Coax Plus Ophthal
5580*	Diagnostix™ Desk Set, PMV Otoscope, Coax Ophthal
55802*	Diagnostix™ Desk Set, PMV Otoscope, Coax Plus Ophthal
5660TN-GAK	Diagnostix™ Global Adapter Kit, EU, UK, AU





USB Rechargeable
Handle Now Avail-
able!



Diagnostic™ 3.5V Portable Instruments

3.5V Diagnostic EENT instruments are among the most critical tools for a comprehensive physical exam.

The Diagnostic™ Portable Set features

- Choose from 3.5V Diagnostic™ Coax or Coax Plus Ophthalmoscope
- Choose from a 3.5V Diagnostic™ standard otoscope or our PMV otoscope
- USB or Plug-In Li-ion Adtronic™ battery handle (see page 8)
- Soft and hard case options depending upon configuration
- 3.5V Throat Illuminator or Dermoscope available
- Choose xenon or AdLED™ lighting technology
- Two-year warranty on instrument heads and handle, lifetime on LED lamp

Component Model	5422 PMV Otoscope	5420 Standard Otoscope	5440 Coax Ophthalmoscope	5442 Coax Plus Ophthalmoscope	5426 Throat Illuminator	5450 Dermoscope	5460 or 5463 Handle	Hard Case	Soft Case
54380* ‡	✓		✓				✓	✓	✓
543802* ‡	✓			✓			✓	✓	✓
54381*	✓						✓		✓
54310* ‡		✓	✓				✓	✓	✓
543102* ‡		✓		✓			✓	✓	✓
54311*		✓					✓		✓
54312*			✓				✓		✓
543122*				✓			✓		✓
54316*					✓		✓		✓
54313						✓	✓		✓

3 Add a 3 to model number to specify USB-C rechargeable handle; leave out when ordering socket plug-in rechargeable handle

* Indicates lamps: X = Xenon, L = LED, A = LED Otoscope & Xenon Ophthalmoscope, B = Xenon Otoscope & LED Ophthalmoscope

‡ Hard case is standard. Add "S" for soft case



Diagnostix™ 2.5V Pocket Instruments

Full-featured pocket sets that sacrifice nothing but size.

The Diagnostix™ Pocket Set features

- Compact portability with a feature-rich design
- AA handle with integrated pocket clip
- Choose from our 2.5V fiber-optic Pocket Otoscope and 2.5V Pocket Ophthalmoscope instrument heads
- Choose xenon illumination for a bright white light and true tissue color or high-performance AdLED™ lighting

- Instruments available individually or in complete two-instrument diagnostic sets
- Available in black, green (GR), royal blue (RB), purple (V), burgundy (BD), or gray (G)
- Offered in either hard fitted or soft nylon storage case
- AA batteries included
- Disposable specula included with Otoscopes
- Spare lamps included with xenon sets (not with LED sets)
- Inspected, assembled, and packaged in the USA from foreign components
- Two-year instrument warranty, lifetime on optics

Ordering Information

Item Description

5110E**	Diagnostix™ Single Handle, Otoscope/Ophthalmoscope
5110N*S**	Diagnostix™ Dual Handle, Otoscope/Ophthalmoscope, Soft Case
5110N* **	Diagnostix™ Dual Handle, Otoscope/Ophthalmoscope, Hard Case
5111N*S**	Diagnostix™ Single Handle, Otoscope, Soft Case
5111N* **	Diagnostix™ Single Handle, Otoscope, Hard Case
5112N*S**	Diagnostix™ Single Handle, Ophthalmoscope, Soft Case
5112N* **	Diagnostix™ Single Handle, Ophthalmoscope, Hard Case

* Indicates lamp Choice: L = LED

** Add BD for Burgundy, G for Gray, GR for Green, RB for Royal Blue, or V for Purple. (Black is standard).

LED Component Model	Xenon Component Model	5120N Otoscope	5140N Ophthalmoscope	5160N Handle	Hard Case	Soft Case 9251XXS
N/A	5110E**	✓	✓	✓	✓	
5110NLS**	5110NS**	✓	✓	✓ (x2)		✓
5110NL**	5110N**	✓	✓	✓ (x2)	✓	
5111NLS**	5111NS**	✓		✓		✓
5111NL**	5111N**	✓		✓	✓	
5112NLS**	5112NS**		✓	✓		✓
5112NL**	5112N**		✓	✓	✓	

** Specify "BD" for Burgundy, "G" for Gray, "GR" for Green, "RB" for Royal Blue, or "V" for Purple when ordering. (Black is standard.)





5211

Proscope™ 2.5V Portable Instruments

ADC's Proscope™ portable instruments offer a wide selection of diagnostic capabilities at budget-friendly pricing.

The ADC 2.5V Diagnostic Set features

- Ophthalmoscope head with lens wheel features 24 corrective lenses from -25 to +40 diopters
- Otoscope head with removable 3X magnifying lens and insufflator fitting
- Bayonet head locking mounts for secure connection to the battery handle
- Welch Allyn disposable polypropylene ear specula in 2.75mm and 4.25mm sizes
- C-cell powered battery handle with rheostatic on/off switch
- Auto-clavable, polypropylene ear specula in 2.5mm, 3.5mm, and 4.5mm sizes
- Adapter to allow use with Welch Allyn or compatible disposable specula
- Sets complete in zippered case with fitted, scratch-resistant, and easy-to-clean insert; single instruments offered in nylon case
- Inspected and packaged in the USA from foreign components
- Batteries NOT included
- One-year warranty



5215



Component Model	5220 Otoscope	5240 Ophthalmoscope	5312-3 Dermoscope	5215-TH Throat Illum.	5211-AD Specula Adapter	5215-NS Nasal Spec.	5215-TD Tongue Dep.	5215-M3 5215-M4 Mirrors	Handle	5211-6 Spec.	Hard Case	Soft Case
5211*	✓				✓				5211-5	✓		✓
5212		✓							5211-5			✓
5312			✓						5260			✓
5210	✓	✓			✓				5211-5	✓	5210-ZC	
5215	✓	✓		✓	✓	✓	✓	✓	5211-5	✓	5215-ZC	

*Also available in veterinary model (5211V) with specula sizes 4, 5, and 7

Popular Adstation Configurations



Shown here: Model 5610-39W
with Advview 2
(Model 9005BPT)

Advview 2 is ordered as
a separate item.

Advview 2

Learn more about the Advview 2 at
www.adctoday.com/advview2

Diagnostix Adstation™ with Advview 2

3.5V Modular Diagnostix Wall System
with Advview 2

- Advview 2 modular diagnostic station with optional pulse oximetry and thermometry
- 3.5V Diagnostix Coax Plus Ophthalmoscope or Coax Ophthalmoscope
- 3.5V Diagnostix PMV Otoscope or Standard Otoscope with specula dispenser
- Add a 429 Non-Contact Thermometer
- Add a 3.5V Diagnostix Throat Illuminator or Dermoscope
- Choose from Halogen/Xenon or high-performance AdLED lighting technology



Shown here:
5680-347
with PMV Otoscope,
Non-Contact Thermometer,
Clock Aneroid
and Wallboard

Diagnostix Adstation™ with Clock Aneroid

3.5V Modular Diagnostix
Wall System with Clock Aneroid

- Diagnostix Clock Aneroid with cuff storage basket, coiled tubing, adult and large adult inflation systems
- 3.5V Diagnostix Coax Plus Ophthalmoscope or Coax Ophthalmoscope 3.5V Diagnostix PMV Otoscope or Standard Otoscope with specula dispenser
- Adtemp 429 Non-Contact Thermometer
- Choose from Halogen/Xenon or high-performance AdLED lighting technology
- Optional wall board for simplified mounting



Shown here:
Model 5680-38W
with e-sphyg 3
(Model 9003-MCC)

**e-sphyg 3 is ordered
as a separate item.**

e-sphyg3™

Learn more about the e-sphyg 3
at www.adctoday.com/esphyg3.

Diagnostix Adstation™ with e-sphyg 3 3.5V Modular Diagnostix Wall System with e-sphyg 3

- E-sphyg 3 NIBP monitor
- 3.5V Diagnostix Coax Plus Ophthalmoscope or Coax Ophthalmoscope
- 3.5V Diagnostix PMV Otoscope or Standard Otoscope with specula dispenser
- Add a 429 Non-Contact Thermometer
- Choose from Halogen/Xenon or high-performance AdLED lighting technology

Diagnostix Adstation™

Basic 3.5V Diagnostix Wall System

- 3.5V Diagnostix Coax Plus Ophthalmoscope or Coax Ophthalmoscope
- 3.5V Diagnostix PMV Otoscope or Standard Otoscope with specula dispenser
- Choose from Halogen/Xenon or high-performance AdLED lighting technology
- Optional wall board for simplified mounting

Can be ordered
as single
instrument.



Shown here:
Model 5613
with Dermoscope



Shown here:
Model 5610-3 with
Standard Otoscope



AtlanticMedicalSolutions.com

704-561-0004

p/n: EENTLIT rev 5
Printed in U.S.A.

