



Hear enhanced high- and low-frequency  
body sounds with **Adjustable Frequency  
Design (AFD)** diaphragms



## Adscope® 608 Convertible Clinician

AFD technology and convertible design make this  
the most versatile of our clinician scopes.

Overall Length: 30"

Standard Weight: 5.8 oz (164 g)\*

### CONVERTIBLE CHESTPIECE

- Convertible chestpiece offers adult diaphragm and bell or adult and pediatric diaphragm configurations
- AFD technology for both adult and pediatric diaphragms; use light pressure for low, bell-type frequency or firm pressure for high, diaphragm-type frequency
- Extra-deep bell for unsurpassed low-frequency response
- Machined from surgical stainless steel\* for outstanding performance and rugged durability

### CLINICIAN HEADSET

- Stainless-steel binaurals\* fixed at a 15° angle to maximize comfort and acoustic seal
- Single lumen construction with reinforcing yoke molded into tubing
- PVC tubing is made without latex or phthalates and stays soft and flexible
- Adsoft™ Plus snap-on silicone eartips (includes both small and large) for exceptional comfort and fit


### ALSO INCLUDES

- Non-chill diaphragm retaining rim for patient comfort
- Accessory kit with additional eartips and scope ID tag

### THE ADC ADVANTAGE

- Assembled, inspected, and packaged in the USA
- Individually acoustic tested in the USA
- Lifetime warranty and free parts forever

\* Weight with pediatric diaphragm is 8.6 oz (243 g). Weight may vary depending on finish. Some special finishes have aluminum binaurals.

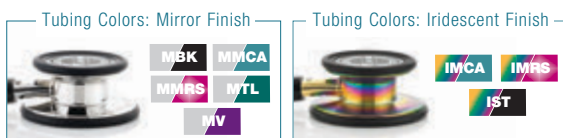

### Ordering Information

Item Description

**608\*\*** Adscope® 608, Convertible Clinician Scope,  
(Specify Color)

**NOTE:** To ensure correct color when ordering,  
add color suffix to model number (e.g., 608**DG**)

