



# Professional and Sensitive Dopplers

## Vascular to Fetal

Newman Medical



ATLANTIC  
MEDICAL  
SOLUTIONS



Over **45 years in the marketplace**; a testament to our longevity and success



**Made in the USA** for superior quality and performance



**5-year all-inclusive warranty** for maximum protection & minimized headache



**Best-in-class customer support**; we always respond quickly

### DigiDop Classic

- Standard size
- Superior battery endurance compared to competitors



### DigiDop II

- Compact size
- Comfortable palm-fitting design
- Equivalent battery endurance to competitors



### DigiDop 901

- Ideal for hospital use
- Clear stereo sound from dual top speakers
- Two probe cradles
- Standard drop-on charging pyramid
- Optional wall mount recharger



#### Vascular Probes



**Model D5**  
5MHz Vascular Probe



**Model D8**  
8MHz Vascular Probe



**Model DPPG**  
Audio TBI PPG Probe



**Model D2**  
2MHz Fetal Probe

#### Fetal Probes



**Model D2W**  
2MHz Waterproof Fetal Probe




**Model D3**  
3MHz Fetal Probe



**Model D3W**  
3MHz Waterproof Fetal Probe

Interchangeable probes can be used on any DigiDop

#### Vascular

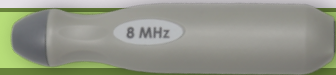


- Great sensitivity
- Ergonomic shape easy to hold, place in sterile glove finger for use in sterile fields
  - **D8 – 8 MHz** - best for shallow vessels like radial, ankle, foot, digits, penile
  - **D5 – 5 MHz** - wider beam, deeper penetration. Good for brachial, carotid, femoral, deep veins

#### Fetal



- Great sensitivity
- Smooth shape, comfortable for patient, easy to hold and clean
  - **D3 – 3 MHz** - wide angle, ideal for early fetal hearts
  - **D3W** - a unique 3 MHz waterproof probe
  - **D2 – 2MHz** - deeper penetration
  - **D2W – 2 MHz** - waterproof

#### DPPG

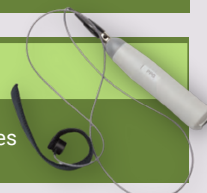

- Unique audio PPG probe for use on digits. Great for toe-brachial studies

#### How to build your Doppler

DOPPLER CODE

PROBE CODE

#### Doppler

DigiDop "Classic"	DD-300	Basic
	DD-301	Basic + Recharge
	DD-700	HR Display
	DD-701	HR Display + Recharge
	DD-901	Tabletop + HR Display + Recharge
DigiDop II	DD-330	Basic
	DD-330R	Basic + Recharge
	DD-770	HR Display
	DD-770R	HR Display + Recharge

#### Probe

D2	2MHz fetal probe
D2W	2MHz waterproof fetal probe
D3	3MHz fetal probe
D3W	3MHz waterproof fetal probe
D5	5MHz vascular probe
D8	8MHz vascular probe
DPPG	Audio PPG probe - includes toe cuff